

## IVORY® Rubber Dam



Made from natural latex, thoroughly washed in the manufacturing stage to reduce surface proteins. Ivory dam has extremely high tear strength, reducing waste from tearing during application. The dam is thin and pliable, allowing easy placement in especially difficult areas. The bright blue or green tint increases visibility by its contrast.

Available in 6x6 or 5x5 (medium weight).



## IVORY® Punch

This precision instrument is made from high quality stainless steel. The unique floating plunger strikes all edges of the cutting disc evenly, guaranteeing a perfect, tear-resistant hole every time. The cutting disc rotates through a range of 6 hole sizes, providing the correct opening for the smallest anterior tooth to the largest molar.

## IVORY® Forceps

The Ivory forceps are made from high grade stainless steel, which will assure many years of use. The forceps will reach around the bow of any clamp without tilting, which is especially important when reversing a clamp. The straight head design allows easy gripping and placement of any clamp size.

## IVORY® Rubber Dam Kits

The Ivory Clamp Starter Kit offers a selection of the most widely used clamps. These popular clamps provide a variety of general use clamps that fulfill most needs. The Rubber Dam Complete Kit is a self-contained unit providing all the necessary equipment for the Rubber Dam Technique. The Complete Kit is especially useful for the practitioner just beginning their practice or for one who is implementing the technique for the first time.

**Kit components include:**

### 8 Clamp Starter Kit

One each of the following clamps:

0 2A 8A 14 2 7 9 14A

Available in Shiny or Matte Finish

### Complete Kit

- (1) Box 6x6 Blue Rubber Dam
- (1) 6x6 Marking Template
- (1) 6 Inch Frame
- (8) Most Popular Clamps (0 2A 8A 14 2 7 9 14A)
- Clamp Forceps
- Rubber Dam Punch

The complete system for tooth isolation.



The complete system  
for tooth isolation.



Ivory®

Please call 1-800-343-5336 for quarterly specials

# Ivory® Rubber Dam Clamps by Heraeus

Heraeus manufactures many styles of uniquely designed IVORY® Rubber Dam clamps to meet the most demanding situations. You will find one of these clamps suitable for practically every application. The tension of each clamp is meticulously hand-set in our factory and is appropriate for the tooth on which it is designed to be used. Ivory Clamps are available in shiny and matte finish as noted in the charts.

Listed below is the complete Ivory Clamp line and suggestions for routine procedures.

## Most Popular Clamps

### PREMOLARS – WINGED

- 00 Small lower, anterior
- 0 Elongated small neck
- 1 General purpose upper
- 2 General purpose lower
- 2A Large neck bicuspid

### PREMOLARS – WINGLESS

- W2 General purpose lower
- W2A Large neck bicuspid

### MOLAR – WINGED

- 3 Small lower
- 4 Small upper for round teeth
- 7 General purpose lower
- 8 General purpose upper
- 8A Upper and lower molar root
- 12A Right partially erupted, 3<sup>rd</sup> molar
- 13A Left partially erupted, 3<sup>rd</sup> molar
- 14 Partially erupted, most popular clamp

### MOLAR – WINGED (CONTINUED)

- 14A Large, partially erupted
- 56 All large molars

### MOLAR – WINGLESS

- W3 Small lowers
- W8A Upper and lower molar root
- W14A Large partially erupted

## Specialty Clamps

### PREMOLARS – WINGED

- 1A Root clamp

### PREMOLARS – WINGLESS

- W0 Elongated small neck
- W1 General purpose upper
- 27N Universal

### MOLAR – WINGED

- 5 Large upper
- 7A Large lower
- 10 Right, small short crown

### MOLAR – WINGED (CONTINUED)

- 11 Left, small short crown
- 26 Small lower molar
- 27 Distal access, large molar

### MOLAR – WINGLESS

- W4 Small upper
- W5 Large upper
- W7 General purpose lower, square teeth
- W8 General purpose upper, round teeth
- W14 Partially erupted

### MOLAR – WINGLESS (CONTINUED)

- W56 Large molar
- 26N Universal upper
- 28N General purpose lower molar

### LABIAL – WINGED

- 6 Central and incisor
- 9 Labial prep, all anteriors

### LABIAL – WINGLESS

- 212 Labial prep, cervical

## Labial

**6**

SHINY



For labial cavities on centrals and cuspids.  
For second bicuspids and second upper molars.

**9<sup>T</sup>**

SHINY  
MATTE



For labial cavities on all anterior teeth.

**212**

SHINY  
MATTE



For labial cavities. Anterior cervical clamp.

**212<sup>SA</sup>**

SHINY




For odd-shaped anteriors and in situations where additional tissue retraction is needed.

# Tooth Diagram

## Tooth Clamp

<b>1</b>	4, 5, 8, 8A, 12A, 14, 26N
<b>2</b>	4, 5, 8, 8A, 12A, 14, 26N
<b>3</b>	4, 5, 8, 8A, 12A, 14, 14A, 26N, 56
<b>4</b>	0, 1, 1A, 2A, 27N
<b>5</b>	0, 1, 1A, 2A, 27N
<b>6</b>	6, 9, 212, 212SA
<b>7</b>	6, 9, 212, 212SA
<b>8</b>	6, 9, 212, 212SA
<b>9</b>	6, 9, 212, 212SA
<b>10</b>	6, 9, 212, 212SA
<b>11</b>	6, 9, 212, 212SA
<b>12</b>	0, 1, 1A, 2A, 27N
<b>13</b>	0, 1, 1A, 2A, 27N
<b>14</b>	4, 5, 8, 8A, 13A, 14, 14A, 26N
<b>15</b>	4, 5, 8, 8A, 13A, 14, 14A, 26N
<b>16</b>	4, 5, 8, 8A, 12A, 14, W14, 14A, 26N, 56

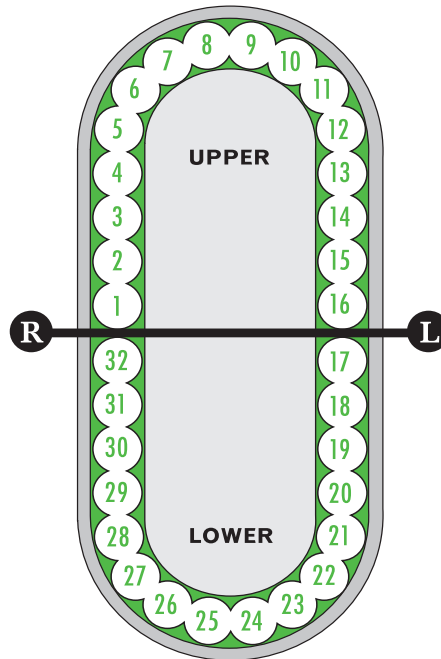

## Tooth Clamp

<b>17</b>	3, 7, 7A, 8A, 12A, 14, 14A, 27, 56
<b>18</b>	3, 7, 7A, 8A, 12A, 14, 14A, 27, 56
<b>19</b>	3, 7, 7A, 8A, 12A, 14, 14A, 27, 56
<b>20</b>	0, 1A, 2, 2A, 6, 27N
<b>21</b>	0, 1A, 2, 2A, 6, 27N
<b>22</b>	00, 6, 9, 212, 212SA
<b>23</b>	00, 6, 9, 212, 212SA
<b>24</b>	00, 6, 9, 212, 212SA
<b>25</b>	00, 6, 9, 212, 212SA
<b>26</b>	00, 6, 9, 212, 212SA
<b>27</b>	00, 6, 9, 212, 212SA
<b>28</b>	0, 1A, 2, 2A, 6, 27N
<b>29</b>	0, 1A, 2, 2A, 6, 27N
<b>30</b>	3, 7, 7A, 8A, 12, 12A, 14, 14A, 27, 56
<b>31</b>	3, 7, 7A, 8A, 12, 12A, 14, 14A, 27, 56
<b>32</b>	3, 7, 7A, 8A, 14, 14A, 27, 56

# Bicuspid

**0°**  
SHINY  
MATTE

For teeth that are elongated or have small necks. High bow, small jaws.  
Also available wingless - W0

**00**  
SHINY  
MATTE

For small lower anteriors and teeth set irregularly in arch. High bow, narrow jaws.

**1°T**  
SHINY  
MATTE

General purpose upper. Jaws fit curvature at gum line.  
Also available wingless - W1

**1A**  
SHINY

Root clamp. Jaws slope downward, grasp roots of bicuspid, anteriors.

**2°T**  
SHINY  
MATTE

General purpose lower. Same as #1, only jaws are flat so as not to impinge on tissue.  
Also available wingless - W2

**2A°T**  
SHINY  
MATTE

For bicuspid with large necks. Broad jaws.  
Also available wingless - W2A

**27N\***  
SHINY

Universal bicuspid clamp. Jaws conform to curvature at gum line.  
Available wingless only

\* Available wingless only  
● Also available wingless  
T Also available with serrated jaws (Tiger clamps)  
A Slightly larger  
N No wing

## Tiger Clamps

1T	Bicuspid, general purpose upper	9T	Anterior labial
2T	Premolar, general purpose lower	14T	Molar, partially erupted
2AT	Premolar, large neck bicuspid	14DT	Molar, distal access, serrated
		56T	For large molars

## Accessories

### IVORY® RUBBER DAM CLAMP STARTER KIT

See page 5 for contents

### IVORY® RUBBER DAM COMPLETE KIT

See page 5 for contents

### IVORY® RUBBER DAM CLAMP BOARD

25 Peg, unnumbered, plastic or stainless steel

### IVORY® RUBBER DAM TEMPLATE 6X6

## Molar

### 3°

SHINY  
MATTE



For small lowers. Narrow, flat jaws so as not to impinge in tissue. Also available wingless - W3

### 4°

SHINY  
MATTE



For small uppers. Narrow jaws designed to fit necks of round-shaped teeth. Also available wingless - W4

### 5°

SHINY



For large uppers. Jaws conform to curvature at gum line. Also available wingless - W5

### 7°

SHINY  
MATTE



General purpose lower. Medium-sized bow and jaws. Flat jaws. For teeth which tend to be square in form. Also available wingless - W7

### 7A

SHINY



For large lowers. Same as #7, only larger. Large bow, broad jaws.

### 8°

SHINY  
MATTE



General purpose upper. For round upper teeth. Jaws fit curvature at gum line, to prevent rotating. Also available wingless - W8

### 8A°

SHINY  
MATTE



For upper, lower molar roots. Four-point jaws grasp root at corners. Jaws set dam beneath gum line when released. Also available wingless - W8A

### 10

SHINY



For small molars with short crowns. Buccal jaw broader than lingual. (10=right; 11=left)

### 11

SHINY



For small molars with short crowns. Buccal jaw broader than lingual. (10=right; 11=left)

### 12A

SHINY  
MATTE



For all lowers, particularly third molars and partially erupted teeth. Available with serrated jaws only. (12A= right; 13A= left)

### 13A

SHINY  
MATTE



For all lowers, particularly third molars and partially erupted teeth. Available with serrated jaws only. (12A= right; 13A= left)

### 14°<sup>T</sup>

SHINY  
MATTE



For partially erupted, undersized, round teeth or for teeth sloping from gum line. Small bow. Also available wingless - W14

### 14A°

SHINY  
MATTE



Same as #14, only larger in every way. Also available wingless - W14A

### 26

SHINY



For small lower molars.

### 26N\*

SHINY  
MATTE



Universal upper with festooned jaws. Available wingless only

### 27

SHINY



For large molar distal cavities.

### 28N\*

SHINY



General purpose lower molar clamp. Available wingless only

### 56°<sup>T</sup>

SHINY  
MATTE



For all large molars. Also available wingless - W56  
56 and W56 are SHINY, 56T is MATTE